

# Wheaton Premium Stoppers

## NEW PRODUCT

# Splash

- Ultra Pure rubber stoppers have low levels of extractables and are latex, nitrosamine, & 2-MBT free
- Complete Coat stoppers are manufactured in a silicone-free environment
- Stoppers available in 2 and 3-Leg Lyophilization, Igloo, and Straight Plug styles

**Ultra Pure Lyophilization Stoppers** are made from a bromobutyl rubber compound containing the highest level of chemical purity. These stoppers demonstrate a very low level of moisture absorption and are ideal for most lyophilization applications. They are compatible with Water for Injection (WFI). They provide a great seal with a low level of extractables.


**Ultra Pure Straight Plug Stoppers** are made from a bromobutyl rubber compound containing the highest level of chemical purity. These stoppers are compatible with Water for Injection (WFI) applications. These stoppers demonstrate ultra low levels of extractables and have a low particulate content. They also remain stable when exposed to gamma irradiation making them very well suited for Ready-for-Sterilization and Ready-for-Use applications.

**Complete Coat Stoppers** are completely coated with a fluorinated polymer that has barrier properties similar to PTFE. The fluorinated polymer coating is the only material that comes in contact with the sample. As a result, chemical reactions due to unwanted interference with the rubber itself are virtually eliminated. These stoppers are manufactured in a silicone-free environment. The specialized manufacturing process produces extremely low levels of visible and subvisible particles. These stoppers provide users the greatest level of confidence in avoiding interference with samples.


All stoppers are packaged 100/shelf pack and 10 shelf packs/case, and they can be purchased from your local laboratory supply distributor.

Cat. No.	For Mouth ID x OD (mm)	Stopper Style	Stopper Material	Shelf Pack	Qty/Case
<b>Ultra Pure Stoppers</b>					
W224100-400	7 x 13	Straight Plug	Gray Bromobutyl 47	100	1000
W224100-401	7 x 13	2-Leg	Gray Bromobutyl 46	100	1000
W224100-402	7 x 13	Igloo	Gray Bromobutyl 46	100	1000
W224100-405	13 x 20	Straight Plug	Gray Bromobutyl 47	100	1000
W224100-406	13 x 20	2-Leg	Gray Bromobutyl 46	100	1000
W224100-407	13 x 20	3-Leg	Gray Bromobutyl 46	100	1000
W224100-408	13 x 20	Igloo	Gray Bromobutyl 46	100	1000
<b>Complete Coat Stoppers</b>					
W224100-420	7 x 13	Straight Plug	Gray Bromobutyl 39	100	1000
W224100-421	13 x 20	Straight Plug	Gray Bromobutyl 39	100	1000

Note: To select the proper size stopper for your serum bottle or vial, use the OD dimension of the stopper, bottle and/or vial.

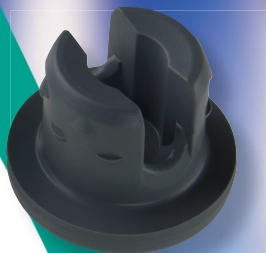

2-Leg Stopper

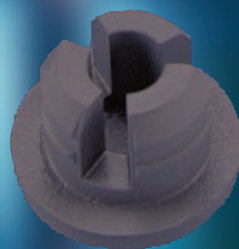

3-Leg Stopper



Straight Plug



Igloo Plug



**WHEATON**

Science Products

a division of WHEATON Industries Inc.

1501 North 10th Street, Millville, NJ 08332-2038  
1.800.225.1437 ■ 1.856.825.1100

[www.wheatonsci.com](http://www.wheatonsci.com)