


WHEATON | SingleScan™ Bar Code Reader

Introducing the **NEW** WHEATON SingleScan Bar Code Reader

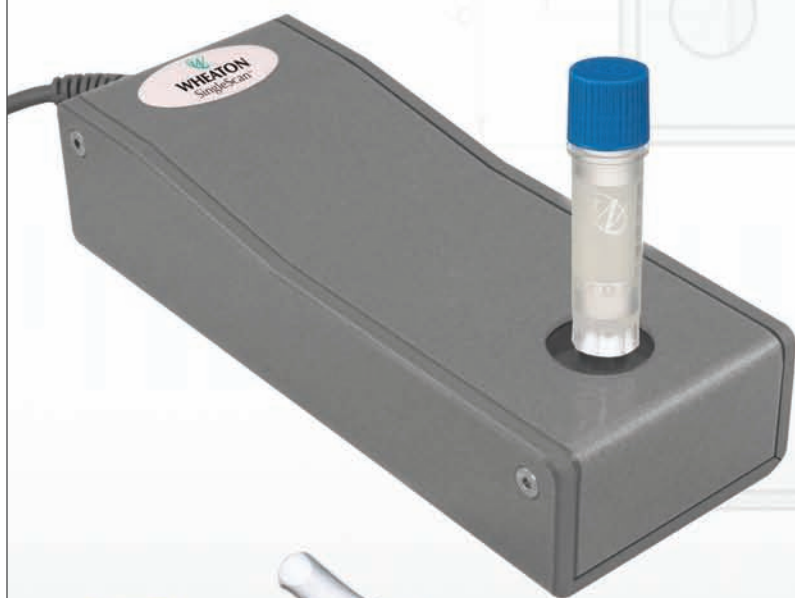

WHEATON | SingleScan Bar Code Reader

The WHEATON SingleScan Bar Code Reader makes sample traceability easy. The SingleScan Bar Code Reader is a sophisticated 2D Data Matrix and 1D Linear Bar Code Reader capable of reading bar codes from a wide range of manufacturers and label printers.

- > Plug and play design enables easy set-up with no software installation required.
- > USB interface allows easy connection to computer.
- > Reads 1D linear bar codes and 2D data matrix bar codes on any vial or ampule.
- > Inputs decoded bar code ID into any software application where cursor is placed.
- > One year warranty.

Cat. No.	Description	Qty
W986000	SingleScan Bar Code Reader	1



> Unique Advantages & Tips

- Bar coding gives you the ability to identify and track your samples, increases the speed of data collection, and eliminate errors which reduces cost.
- It is critical that bar codes are read correctly every time, remain adhered to the container, and are durable enough to withstand severe environments.
- WHEATON can customize a bar code solution to meet your specific application, while guaranteeing no duplicate numbers and sequence integrity.

See Related Accessories on Reverse Side

Related Accessories

CryoELITE™ Cryogenic Vials

- > CryoELITE Cryogenic Vials are the perfect solution for storing, tracking and transporting your irreplaceable samples.
- > Unrivaled cap seal exceeds DOT and IATA regulations ensuring ultimate protection of samples during transportation and demanding freeze-thaw handling.
- > Loctagon™ vial skirt provides stability in the freestanding position; locks into the CryoELITE Benchmate Rack to provide easy open and close with one hand.
- > E-Z Ex-Traction® conical well design allows precise sample retrieval for both automated and manual systems.
- > 7 cap colors for convenient sample identification.
- > Lot certified RNase/DNase and endotoxin free, sterile and non-pyrogenic.




CryoELITE Bar Code Inserts

- > Optional 2D Data Matrix Bar Code Insert allows for immediate bar coding of your samples.
- > Inserts can also be applied to the vial at a future date, which eliminates jeopardizing the integrity of your sample by transferring it to another vial.
- > The bar code insert press fits and locks into in the bottom of the CryoELITE freestanding vial.
- > Bar code number sequencing guaranteed to not duplicate.



KeepIT™ Freezer Boxes

- > KeepIT 100 accommodates up to 100 internal threaded (1.2 - 2mL) cryogenic vials.
- > KeepIT 81 accommodates up to 81 external threaded (1.2 - 2mL) cryogenic vials.
- > Convenience of sample storage and scanning capability in one box.
- > Made from Eastman Tritan™ BPA free shatter resistant resin.
- > Color coded to match CryoELITE vial cap colors for fast visual batch identification.
- > Holes in bottom facilitate scanning 2D Data Matrix Bar Code Inserts.

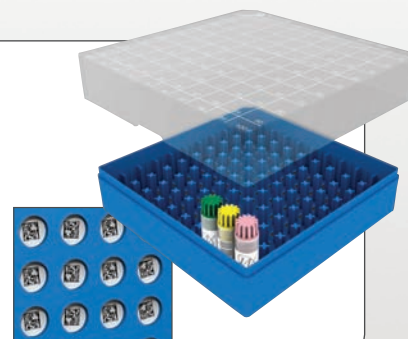

Custom Bar Coding

- > Customized bar code solutions to meet your specific application.
- > Choose from a wide range of adhesive, etching, and print technologies to ensure accurate identification and tracking.
- > 100% scannable.
- > Packaged and sequentially numbered to your specifications.

